



# 1602 Essentia

Surgical Table



# Providing Essential Needs at a High Value



The 1602 Essentia provides exclusive top rotation in a streamlined design. Achieve surgical and imaging access with its 180-degree rotation similar to our specialty tables at an impressive value. The 1602 Essentia is designed for performance and helps any operating room or surgery center manage more procedures with one economical table.

## Exceptional Imaging Efficiency

- Achieve extensive imaging access at the head and foot end of the table
- Easily permit rotational access by a turn of the handle
- Simple rotation toward the C-arm station at the foot end of the table

## Simple Surgical Positioning

- Great capacity to handle increased specialty cases
- Seamless integration with accessories and imaging equipment
- Ready to perform for a wide range of patients and procedures

## Reliability and Versatility

- Electro hydraulic cylinders for quiet and smooth operation
- Battery powered backup in case of power loss
- One-touch Pendant Control for simple maneuvering



Value-Driven  
Design



Versatile  
Positioning  
Options



Increased  
Patient Safety



## 1602 Essentia Series



### 1602 Essentia Details:

- 180° top rotation
- 500 lb lift and articulation capacity
- Battery power
- Positioning options include 90° back-up position
- Durable stainless steel base is easy to clean
- 27 <sup>3</sup>/<sub>16</sub>" - 40 <sup>15</sup>/<sub>16</sub>" adjustable table height
- 23° lateral tilt
- 28° Trendelenburg and reverse Trendelenburg
- 90° back-up for beach chair and fowler positioning

# Performance & Value for Versatile Operating Rooms

The 1602 Essentia general purpose surgical table works to fit your unique procedures covering essential needs like:

- Bariatric
- Cardiovascular
- ENT/ Ophthalmology
- Imaging
- Robotic
- Urology
- and more!



Self-leveling brakes for stability on uneven floors



Optional 5" side extensions widen the table surface to 30"



Power cord can be easily removed when using battery power



27 <sup>3</sup>/<sub>16</sub>" Minimum Table Height



40 <sup>15</sup>/<sub>16</sub>" Maximum Table Height



30° Flex



90° Reflex





0° Leg Up



100° Leg Down

# Simple, One-Touch Pendant Control

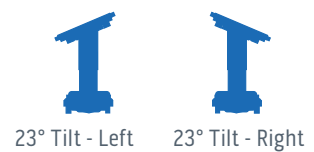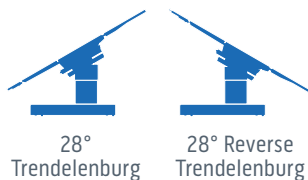

# Positioning Guide



Upper Body Imaging



Lower Body Imaging



Spinal



Vascular/Endovascular



Ophthalmic/ENT



Kidney/Thoracic

# Positioning Guide



Neuro (Neck)



Neuro (Lumbar)



Shoulder Arthroscopy



Cysto/Gyn



Gall Bladder/Thoracic



Bariatric

# Skytron's Solutions



## Sterile Processing

Processing Sinks and Stainless Offerings  
Washers, Ultrasonics, and Cart Washers  
Prep & Pack Tables & Sterilizers



## Clinical

Surgical Tables  
Booms & Lights  
Full Offering of Stainless Products



## Architectural & Environmental

Modular Walls  
Configurable Window and Door Offerings  
Built-in Cabinetry



## Clinical Business Intelligence

Video Integration  
Asset Tracking & Supply Chain Solutions  
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

**To learn more, visit our website at [www.skytron.com](http://www.skytron.com)**



5085 Corporate Exchange Blvd. SE  
Grand Rapids, MI 49512  
1.800.SKYTRON (759-8766)  
email: [info@skytron.com](mailto:info@skytron.com)  
[www.skytron.com](http://www.skytron.com)



REV3 03-24  
PN L2-010-109