



3503 EZ Slide

Surgical Table

The Standard in Full Body Imaging & Patient Positioning



Imaging and personalized height adjustment have never been easier with the 3503 EZ Slide. Its market-leading 21" (53 cm) of top slide, allows for more movement of the patient through the imaging window. Also, attain comfortable working access for a broad range of specialties with 3503's height range up to 43" (109 cm). The 3503 EZ Slide allows you to do more with fewer constraints with a lower cost of ownership.

Increased Imaging Efficiency

With 21" (53 cm) of top slide, the 3503 EZ Slide makes full body imaging convenient. A simple click of a button effortlessly moves the patient toward the C-arm, eliminating the need for risky repositioning. The table also features a slim profile base that stays clear of mobile imaging equipment.

Surgical Positioning Flexibility

- Removable leg section for simple setup and positioning
- Quickly integrates with accessories like the 40" (102 cm) carbon fiber extension
- Made with sturdy, standardized side rails acceptable to most accessories
- Preprogrammed positioning options created by the touch of a button

Improved Patient Safety

- Achieve top slide at full weight capacity at 600lbs (272kg)
- No confusion on weight limits in articulation



Comfortable
Surgical
Access



Powerful
Articulation
Capacity



One-Touch
Imaging
Access



3503 EZ Slide Series

21" (53 cm) Top Slide, General Purpose Table

Providing 21" (53 cm) of top slide, the 3503 EZ Slide takes the work out of imaging setup. At the click of a button, the patient is moved through the expansive imaging window.



Differentiating Features


- 21" (53cm) top slide
- 7½" (19cm) slide to head and 13" (33 cm) slide to foot
- 700 lb (318kg) lift and 600 lb (272kg) articulation weight capacity
- Removable head and leg sections for quick accessory attachment and positioning
- Durable stainless steel base is easy to clean
- 25" - 43" (64 cm - 109 cm) adjustable table height
- 20° lateral tilt
- 30° Trendelenburg and 30° reverse Trendelenburg



Convenient and Comfortable Imaging Access

The 3503 EZ Slide will offer 21" (53 cm) of top slide at its full articulation weight capacity. Designed to give comfortable, waist-height access to taller surgeons, the table offers our maximum table height of 43" (109 cm). As the table top slides toward the foot end, the head section can be conveniently detached and repositioned as a foot extension.

The 3503 EZ Slide also brings confidence to the operating room. Recognized for their reliable performance, Skytron surgical tables are consistently rated as having a lower cost of ownership.

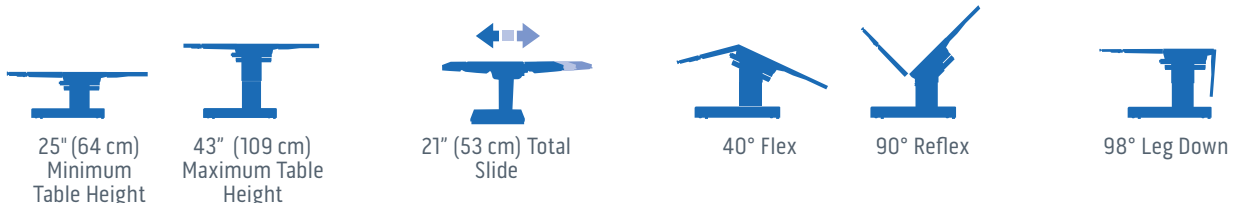



Back lit one-touch positioning functions include Trendelenburg, tilt, and more




Battery life/diagnostic indicator



Easy to use clip style power cord



Back-up controls for emergency situations

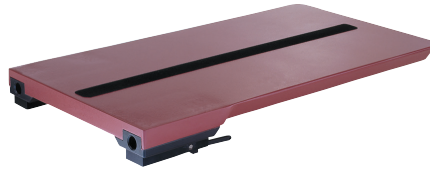

Optional 40" (102 cm) carbon fiber leg section for lower body imaging



Self-leveling brakes for stability on uneven floors

One-Touch Pendant Control

- Simple one-touch operation
- Backlit control for high visibility in low lighting
- Audible alarm for kidney bridge, back up articulation and top slide obstruction
- With PSR (Patient Safe Return)
- Recessed buttons prevent accidental activation
- Durable and water resistant
- Can be used in conjunction with optional three function foot control



90° Back Up



40° Back Down



30° Trendelenburg



30° Reverse Trendelenburg



20° Tilt - Left



20° Tilt - Right

Positioning Guide



Upper Body Imaging



Lower Body Imaging



Lap Nissen



Spinal



Ophthalmic/ENT



Nephrectomy

Positioning Guide



Neuro (Neck)



Neuro (lumbar)



Shoulder Arthroscopy



Urology



Lower Extremity



Bariatric

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Configurable Window and Door Offerings
Built-in Cabinetry



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV2 03-24
PN L2-010-140