



6302 Elite

Surgical Table

Bringing value and efficiency to every OR.



Offering exclusive top rotation with an efficient design, the 6302 Elite will bring enduring value to any operating room. It's 210° table top rotation offers simple and quick setup and accommodates imaging requirements with ease. With efficient operation and market-leading reliability, the 6302 Elite will continue to deliver unsurpassed value long after the initial purchase.

Imaging Efficiency

With the simple turn of a handle, the 6302 Elite can rotate up to 210° – making imaging access easy at the head and foot end of the table. The table provides hassle-free, full-body imaging with its exclusive rotation that allows effortless setup with a C-arm. Without any complex setup requirements, any operating room can reduce the time spent preparing for imaging cases.

Surgical Positioning Flexibility

Make every procedure as efficient as possible with 210° top rotation. The patient can seamlessly transition between the anesthesiologist and surgeon during head-end procedures. With a power beach chair positioner and a removable leg section for quick accessory attachment, you can save valuable time.

Improved Patient Safety

Because the 6302 Elite table top turns independently of the base, the table brakes can remain firmly in place throughout rotation making patient maneuvers safe and quick. In addition, the table's impressive track record of reliability gives you peace of mind that it will continue to deliver unmatched results.



Impressive
Value



Versatile
Positioning
Options



Efficient
Top
Rotation



6302 Elite Series

210° Top Rotation, General Purpose Table

Skytron offers true top rotation allowing the surgeon or C-arm to remain at the foot end of the table where the base does not interfere. With a removable leg section, accessory attachment is easy. Use with our optional 40" carbon fiber extension for improved imaging and quick setup.

Differentiating Features

- 210° top rotation
- 1,000 lb lift and 600 lb articulation weight capacity
- Removable leg section for quick accessory attachment and positioning
- Power beach chair positioner for one-touch set-up
- Durable stainless steel base is easy to clean
- 25"-39 7/8" adjustable table height
- 30° lateral tilt
- 30° Trendelenburg and reverse Trendelenburg





Value-Driven Efficiency for Every Procedure

Affordable and Dependable

- Great table for the value-driven consumer - unique 210° rotation feature provides better imaging and surgical access safer than reversing the patient
- Very low down time - you can depend on the durability and sustainability of our surgical tables

ENT

- No table base under patient - easier access
- Larger area under patient - perform more surgeries while seated

Imaging

- Top rotation - allows more space for C-arm
- Removable leg section – allows for Carbon fiber extension accessory - completely radiolucent from side to side

Bariatric

- 1,000lb lift - exceeds most bariatric needs
- 600lb articulating capacity - easily obtain patient position

Urology

- 23" of imaging length - simplifies imaging of the kidneys, ureter, and bladder
- Lightweight transfer board accessories - easy for staff to maneuver



25" Minimum
Table Height



39 7/8" Maximum
Table Height



40° Flex



90° Reflex



20° Leg Up



100° Leg Down



One-touch positioning functions include return-to-level and beach chair



Battery life/diagnostic indicator



Easy to use clip style power cord



Back-up controls for emergency situations



Optional 40" carbon fiber leg section for upper and lower body imaging



Self-leveling brakes for stability on uneven floors

One-Touch Pendant Control

- Simple one-touch operation
- Backlit control for high visibility in low lighting
- Anatomically correct return-to-level
- Anatomically correct beach chair position
- Recessed buttons prevent accidental activation
- Durable and water resistant
- Can be used in conjunction with optional three function foot control



90° Back Up



40° Back Down



30° Trendelenburg



30° Reverse Trendelenburg



30° Tilt - Left



30° Tilt - Right

Positioning Guide



Upper Body Imaging



Lower Body Imaging



Lap Nissen



Urology



Ophthalmic/ENT



Nephrectomy




Neuro (Neck)



Neuro (Lumbar)



Shoulder Arthroscopy

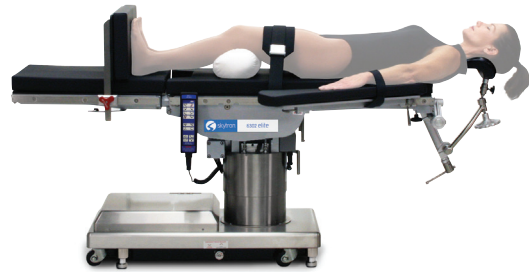

Neuro (Head)



Hip Revision



Cardiovascular

For a full offering of Table Accessories please visit skytron.com/surgical-table-accessories

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Configurable Window and Door Offerings
Built-in Cabinetry



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service. Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions. **To learn more, visit our website at www.skytron.com**



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV3 03-24
PN L2-010-80