



6702 Hercules

Surgical Table



Giving operating rooms the ability to do more.



Combining impressive 1,200 pound lift capacity with exclusive table top rotation, the 6702 Hercules exceeds all standards in general purpose surgical tables. Skytron offers 210° table top rotation – making imaging and patient positioning easy and safe. From pediatric to bariatric, the 6702 Hercules gives any operating room the versatility to handle more procedures and specialties.

Imaging Efficiency

At the turn of a handle, the 6702 Hercules allows superior imaging access at both the head and foot end of the table. The table can quickly rotate toward the C-arm and the low-profile base stays clear of the imaging equipment. Skytron takes the work out of imaging setup.

Surgical Positioning Flexibility

A remarkable 1,200 pound lift and 1,000 pound articulation capacity makes the 6702 Hercules ready for patients of diverse sizes. A low 23⁵/₈" table height provides great access to larger patients. The removable back and leg sections allow seamless integration with accessories such as the power beach chair positioner – providing automatic, one-touch setup for shoulder surgery.

Improved Patient Safety

With a table top that rotates independent of the base, the 6702 Hercules allows the brakes to remain firmly in place while moving the patient and maximizing surgical access. This unique feature sets the 6702 Hercules apart and makes any operating room safe.



Increased
Surgical
Specialties



Generous
Weight
Capacity



Efficient
Top
Rotation



6702 Hercules Series

210° Top Rotation, General Purpose Table

Skytron offers true top rotation, allowing the surgeon or C-arm to remain at the foot end of the table where the base does not interfere. And with the removable back and leg section, the power beach chair can be attached for the most comfortable shoulder surgery setup available.



Differentiating Features

- 210° top rotation
- 1,200 lb lift and 1,000 lb articulation weight capacity
- Removable back and leg sections for quick accessory attachment and positioning
- Power beach chair positioner for one-touch set-up
- Durable stainless steel base is easy to clean
- 23 5/8"-41 1/2" adjustable table height
- 30° lateral tilt
- 30° Trendelenburg and reverse Trendelenburg



One Table, More Specialties

Orthopedics – Shoulders

- 180° top rotation - position closer to surgical site
- Removable back section - allows attachment for table powered beach chair

ENT

- No table base under patient - easier access
- Larger area under patient - perform more surgeries while seated

Imaging

- Top rotation - allows more space for C-arm
- Carbon fiber extension accessory - completely radiolucent from side to side

Bariatric

- 1,200lb lift - exceeds most bariatric needs
- 1,000lb articulating capacity - easily obtain patient position

Urology

- 23" of imaging length - simplifies imaging of the kidneys, ureter and bladder
- Light weight transfer board accessories - easy for staff to maneuver



23⁵/₈" Minimum
Table Height



41¹/₂" Maximum
Table Height



40° Flex



90° Reflex



20° Leg Up



100° Leg Down



One-touch positioning functions include return-to-level and beach chair



Battery life/diagnostic indicator



Easy to use clip style power cord



Optional side extensions widen the table top to 30"



Optional 40" carbon fiber leg section for upper and lower body imaging



Self-leveling brakes for stability on uneven floors

One-Touch Pendant Control

- Simple one-touch operation
- Backlit control for high visibility in low lighting
- Anatomically correct return-to-level
- Anatomically correct beach chair position
- Recessed buttons prevent accidental activation
- Durable and water resistant
- Can be used in conjunction with optional three function foot control



90° Back Up



40° Back Down



30° Trendelenburg



30° Reverse Trendelenburg



30° Tilt - Left



30° Tilt - Right

Positioning Guide



Upper Body Imaging



Lower Body Imaging



Lap Nissen



Urology



Ophthalmic/ENT



Nephrectomy

Bariatric Options



Abdominal



Gall Bladder




Neuro (Neck)



Neuro (Lumbar)



Shoulder Arthroscopy

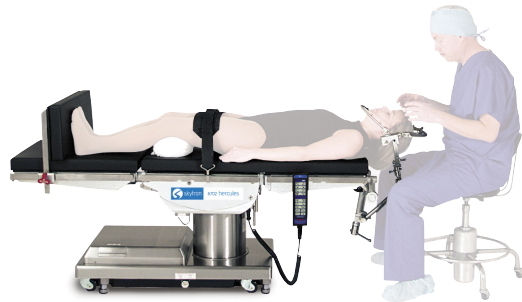

Neuro (Head)



Hip Revision



Cardiovascular



Cysto/Gyn

For a full offering of Table Accessories please visit skytron.us/products/pages/table_accessories

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Configurable Window and Door Offerings
Built-in Cabinetry



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service. Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.
To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV3 03-24
PN L2-010-77