



GS70 Salus

Guardian Series Surgical Table



Setting the Standard for Safety in the OR

The GS70 Salus premier surgical table — the next level of surgical tables from Skytron. Designed with patient and staff safety in mind, the GS70 offers unique and advanced features which will keep any OR informed and safe. Named for Salus, the Roman goddess of safety and well-being, the GS70 Salus is the first of the Guardian Series surgical tables from Skytron. Be more informed and aware before, during, and after surgery. Skytron delivers a low cost of ownership and provides exceptional support for the life of the product.




Verbal Alerts
Increase Patient
Safety



Intelligent
Accessories



Lower Risk
with Data
Capture

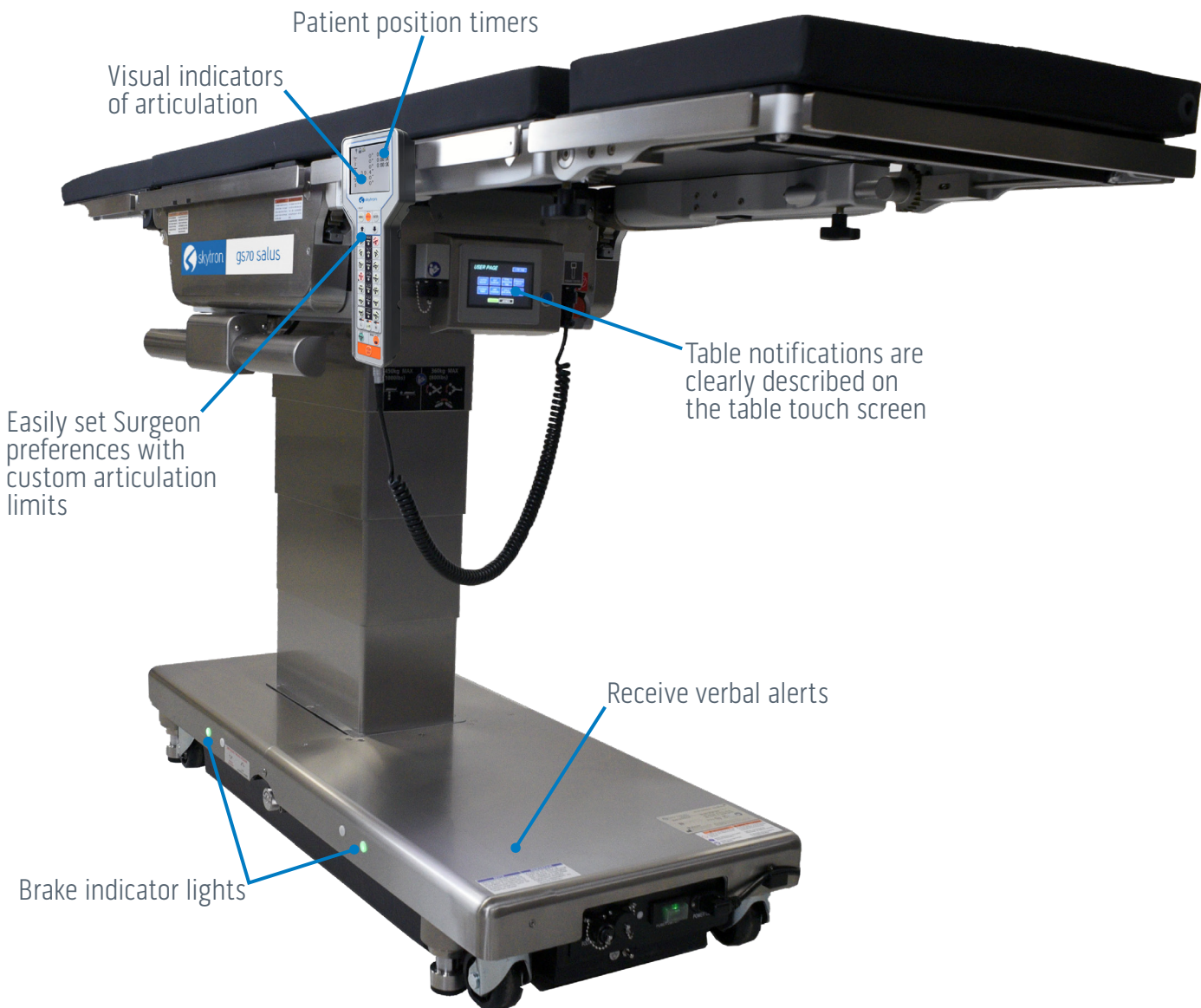

Communication is Everything

High performance teams communicate well. Operating Rooms are uniquely high-risk and time-sensitive locations, with numerous distractions, barriers, and challenges. The GS70 Salus surgical table is designed to support great communication among high-performing teams.

- Easily set Surgeon preferences with custom articulation limits
- Improve patient safety with a verbal alert when positioning may create a tipping hazard
- Increase safety with intelligent accessories that reduce collisions
- Achieve perfect patient positioning with visual articulation graphics (Trendelenburg, Tilt)
- Support best practices and continuous improvement with patient position timers
- Immediate warning for table notifications are clearly described on the table touch screen
- Support research and continuous improvement using captured performance data



Performance Counts

Skytron's general surgical tables add value by supporting a broad range of cases, enabling the most flexibility in using precious Operating Room resources.

Smart accessory attachment

- Automatically adjusts articulation parameters
- Recognizes select accessory engagement

Surgical and Imaging Access

Industry leading table slide and articulation

Exceptional patient positioning options


High articulating weight limits with industry leading top slide

Future proof

Tremendous flexibility and versatility to support any general OR procedure


- Powered split-leg
- Removable head, back, and leg sections
- One touch beach chair articulation
- Wired and wireless pendant control
- Full offering of accessories

One-Touch Pendant Control

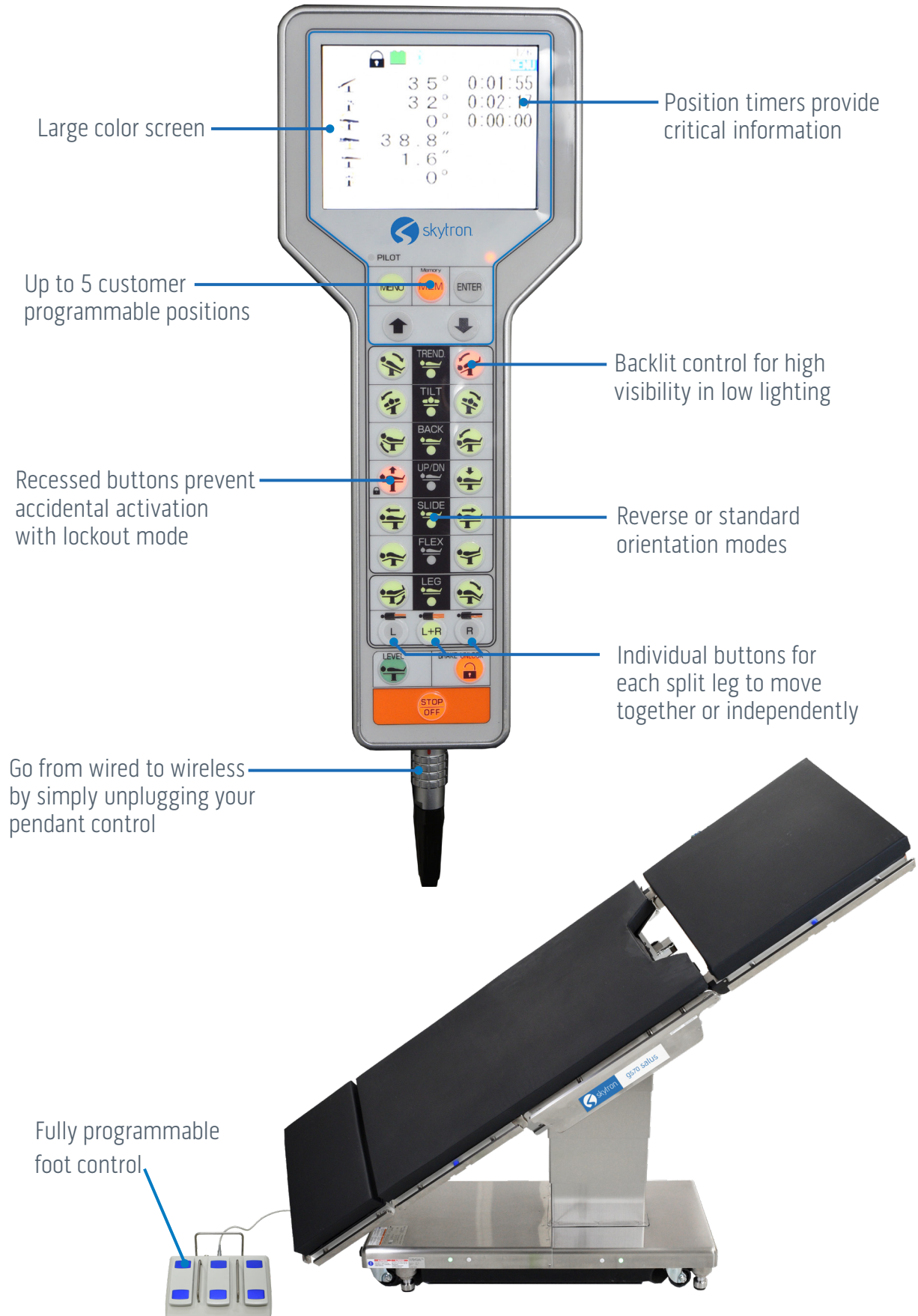





Updated release levers increase safety and allow full table slide in any position.

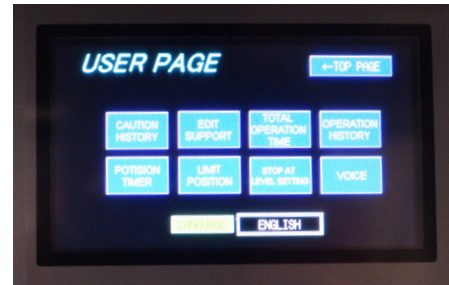

Easy to read, high-contrast, color touch screen for backup control and over 40 alert codes that can be visualized for service diagnostics.

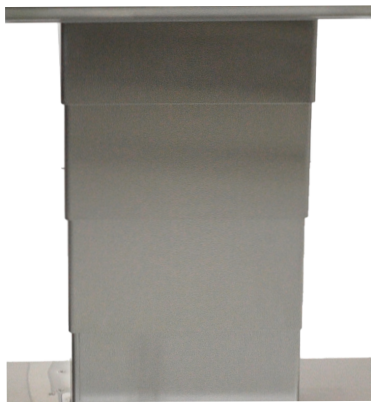

Skytron elevation columns set the standard for a solid, stable platform for Surgeons.

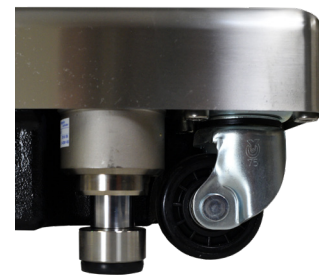

Weight sensing brakes provide stability information to keep your patient and staff safe.

Quality Matters

Skytron clinically superior surgical tables are built from the ground up with long-term, safe operation in mind.

- Best brakes in the business – designed for exceptional strength and durability
- Brake lock indicators – giving teams clear visual checks to see if the table is safe for the patient and staff
- Thick, strong stainless steel lasts for years without cracks and stands up to tough cleaning
- Flush fasteners limit catching of table drapes or pads, simplifying table preparation
- Skytron elevation columns set the standard for a solid, stable platform for Surgeons, and a safe environment for patients
- Hydraulic valves designed specifically for surgical table articulation – not from the hardware store
- Less downtime – a table tagged for non-use and waiting for parts is of little value
- Total cost of ownership – Skytron is the leader in 10 year total cost of ownership

Differentiating Features

- Program up to five user-defined positions
- New accessory connections ensure a solid and secure connection between the accessory and the table
- 26" (66 cm) top slide
- 10" (25 cm) slide toward head and 16 ½" (42 cm) slide toward foot
- 1,000 lb (454kg) lift and 800 lb (363kg) articulation weight capacity
- One-touch setting for beach chair and flex
- Durable stainless steel base is easy to clean
- 35° lateral tilt
- 40° Trendelenburg and 30° reverse Trendelenburg
- Adjustable table height
 GS70 20 ½" to 39 ½" (52 cm to 100 cm)
 GS70H 23" to 41 ½" (58 cm to 105 cm)



Numbers That Lead the Industry

800 Pounds
(363kg)
Articulation Weight Capacity

1000 Pounds
(454kg)
Lift Weight Capacity

Durable 304
Stainless Steel

26" (66 cm)
Top Slide

GS70
20 ½" (52 cm)
Minimum
Table Height

GS70
39 ½" (100 cm)
Maximum Table
Height

GS70H
23" (58 cm)
Minimum
Table Height

GS70H
41 ½" (105 cm)
Maximum Table
Height

40° Flex

90° Reflex

20° Leg Up

100° Leg Down

90° Back Up

40° Back Down

40°
Trendelenburg

30° Reverse
Trendelenburg

35° Tilt - Left

35° Tilt - Right

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Configurable Window and Door Offerings
Built-in Cabinetry



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV1 03-24
PN L2-010-160